
cecotec

NANA BABYCARE SMARTWEIGHT

Báscula de bebé./Baby scale

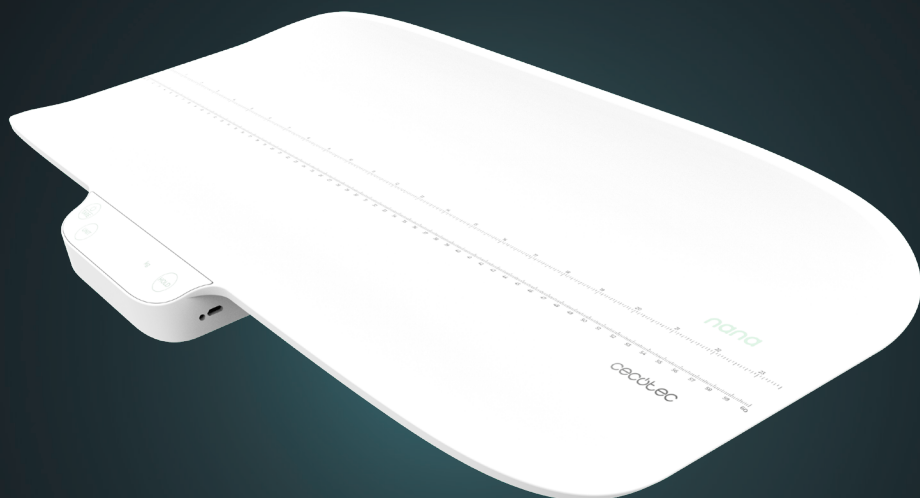Manual de instrucciones
Instruction manual
Manuel d'instructions
Bedienungsanleitung
Manuale di istruzioni
Manual de instruções
Handleiding
Instrukcja obsługi
Návod k použití

Instrucciones de seguridad	4
Safety instructions	6
Instructions de sécurité	8
Sicherheitshinweise	10
Istruzioni di sicurezza	12
Instruções de segurança	14
Veiligheidsinstructies	16
Instrukcje bezpieczeństwa	18
Bezpečnostní pokyny	20

ÍNDICE

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SAFETY INSTRUCTIONS

Read these instructions thoroughly before using the device. Keep this instruction manual for future reference or new users.

- This appliance is intended for household use and not for commercial purposes.
- Keep the device and its packaging materials out of the reach of children and unsupervised pets.
- Do not use the appliance in case of visible damage. In case of malfunction or breakdown, please contact the official Cecotec Technical Support Service.
- Keep the appliance away from strong heat sources.
- Do not use the appliance outdoors.
- Do not expose the device to extreme temperatures, moist environments, or direct sunlight.
- Wide temperature variations could cause condensation. Should this be the case, wait until the appliance is at room temperature before switching it back on.
- Do not drop the device or jump on its measuring surface.
- Supervise children to prevent them from playing with the appliance.

Instructions on batteries

- Do not burn or expose the batteries to high temperatures, as they may explode.
- Both batteries and cells can leak under extreme conditions. Should they leak, do not touch the liquid. If the liquid gets into contact with skin, wash immediately with soap and water. If the liquid gets into the eyes, wash them immediately with clean water for a minimum of 10 minutes and seek medical attention. Wear gloves when handling the

battery or cells and dispose of it immediately in accordance with local regulations.

- Avoid contact between batteries/cells and small metal objects such as paper clips, coins, keys, nails, or screws.
- Do not tamper with batteries or cells in any way.
- Batteries must be removed for safe disposal.
- Keep batteries out of the reach of children when there is a risk of ingestion.
- Battery ingestion can cause burns, soft tissue perforation, and death. It can cause severe burns within two hours of the ingestion.
- In case of battery ingestion, please seek medical help immediately.
- The batteries' power supply terminals and the cells must not be short-circuited.
- If the appliance is to be stored for a long period of time, it is advisable to remove the batteries.
- Batteries and cells must be inserted in their correct position. When replacing batteries, check the polarity.

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1. PARTS AND COMPONENTS

Fig. 1

1. Weighing surface
2. Buttons
3. Display
4. Measurement scale

NOTE:

The graphics in this manual are schematic representations and may not exactly match the device.

2. BEFORE USE

- This appliance is packaged in a way as to protect it during transport. Take the device out of its box and remove all packaging materials. You can keep the original box and other packaging elements in a safe place. This will help you prevent damage to the device when transporting it in the future. In case the original packaging is disposed of, make sure all packaging materials are recycled accordingly.
- Make sure all parts and components are included and in good conditions. If there is any piece missing or in bad conditions, contact the official Cecotec Technical Support Service immediately.

Box content

- Scale
- Instruction manual
- Batteries 2 x AAA

3. OPERATION

When using a baby scale, it is important to follow some guidelines and precautions to ensure accurate measurements and the safety of the baby. Here are some suggestions:

- Place the scale on a firm and level surface. Avoid placing it on carpets or soft surfaces that may affect the accuracy of measurements.
- Avoid damaging, dropping, or knocking the product, as this may affect its accuracy.
- Do not place or use the product in damp, very hot or very cold places.
- This scale is for domestic use only, do not use it for professional purposes.

- Place the scale in a place where there are no drafts that may affect the stability of the scale or the comfort of the baby.
- Make sure there is enough space around the scale so that you can place and remove the baby without difficulty.
- Keep the scale in an environment with a suitable temperature to avoid possible fluctuations in measurements.
- Never leave your baby alone on the scale. Always maintain constant supervision while the baby is being weighed.
- Do not hold or force the baby while on the scale. Allow the baby to be in a comfortable and natural position.
- Make sure you know the manufacturer's recommended weight limit. Do not exceed this limit to avoid damage to the scale and to ensure accurate measurements.
- Use the tare function according to the manufacturer's instructions to obtain accurate measurements by subtracting the weight of clothing or other objects.
- The scale is battery operated, make sure it has sufficient power for operation.
- Before placing the baby on the scale, check its stability and make sure that no parts are loose or damaged.

By following these guidelines and precautions, you can use the baby scales safely and get accurate measurements of your baby's weight. It is always advisable to read and follow the specific instructions provided.

Buttons

Fig. 2

1. ON/OFF - TARE button

On/Off function

Press and hold this button to switch the scale on or off as desired.

TARE function

The tare function allows the scale to automatically subtract the weight of the nappy, clothing or other objects that may be present when weighing the baby, thus only displaying the net weight of the baby.

First, place the accessory on the scale and press the ON/OFF - TARE button. The scale will set the current weight as zero (Fig. 3). Then, place the baby next to the accessory, the scale will only show the difference in weight, i.e. the weight of the baby. This will make it easier to monitor the baby's growth and development more accurately and comfortably.

2. HOLD button

Press this button to maintain the weight shown on the display.

ENGLISH

When the weight displayed on the scale remains stable, press the HOLD button so that the figure does not disappear from the display. If you want to cancel this function, press the button again.

With the HOLD function active, the weight will remain on the display, even if you remove the baby from the scale.

3. UNIT button

Press this button repeatedly before weighing to select the desired weight unit on the display. The available weight units are g-kg-lb-oz-lb: oz (Fig. 4).

Use of the scale

(Fig. 5)

1. Place the scale on a table, desk, or other flat, hard surface. Do not place it on carpets.
2. Press the ON/OFF button to switch it on. Press the UNIT button to change the desired unit of weight. Next, press TARE and wait until "0000" appears on the display.
3. Place the baby on the scale. The baby's weight will appear on the display in the weight unit you have selected.
4. Press the HOLD button if you want the weight to stay on the display and not disappear. If you want to cancel this function, press this button again.

Automatic switch-off function (battery saving)

The scale will automatically switch off if it has not been used for two minutes.

Placing the batteries

1. Open the battery compartment at the back of the scale.
NOTE: please note the polarity of the batteries before inserting them (Fig. 6).
2. Position them correctly according to the instructions.

4. CLEANING AND MAINTENANCE

- Never immerse the appliance in water or other liquids. No liquid must enter the appliance during cleaning.
- Clean the appliance with a slightly dampened, soft cloth and, if necessary, a mild cleaning agent.
- Do not use any caustic or abrasive cleaning agents to clean the device, as these can damage the surface.

WARNING:

- Remove the batteries if the scale is not to be used for a long period of time.
- When not in use, store the appliance in a dry, clean place which is protected from direct sunlight and is out of the reach of children.

5. TROUBLESHOOTING**Error codes**

Replace the batteries when "LO" appears on the display.

If the message "ERR" appears on the display, it means that the scale is overloaded (the maximum weight it can support is 40 kg). Remove weight from the scale.

Please note that the minimum weight is 50 g, the scale will not detect lower weights and will display "0000".

Troubleshooting table

Problem	Possible cause	Possible solution
The scale does not switch on	Batteries exhausted or incorrectly fitted.	Check the installation of the batteries and replace them if necessary.
Inaccurate measurements	Unlevel or uneven surface.	Place the scale on a flat and level surface. Check that there are no objects under the scale.
Inconsistent results	Electromagnetic interference.	Move the scale away from electronic devices that may cause interference.
Tare error	Extra object on the measuring surface.	Remove any additional objects before weighing the baby.
Fluctuating measurements	Unstable environmental conditions.	Place the scale in a place with a stable temperature and no draughts.

6. TECHNICAL SPECIFICATIONS

Product reference: 01902

Product: Nana BabyCare SmartWeight

Minimum weight: 50 g

Maximum weight: 40 kg

Power supply: 2 x 1.5 V AAA batteries.

Range of error:

Up to 5 kg: ± 30 g

From 5 kg to 20 kg: ± 50 g

From 20 kg to 40 kg: ± 100 g

Technical specifications may change without prior notification to improve product quality.

Made in China | Designed in Spain

7. DISPOSAL OF OLD ELECTRICAL AND ELECTRONIC APPLIANCES

This symbol indicates that, according to the applicable regulations, the product and/or batteries must be disposed of separately from household waste. When this product reaches the end of its shelf life, you should dispose of the cells/batteries/accumulators and take them to a collection point designated by the local authorities.

Consumers must contact their local authorities or retailer for information concerning the correct disposal of old appliances and/or their batteries.

Compliance with the above guidelines will help protecting the environment.

8. TECHNICAL SUPPORT AND WARRANTY

Cecotec shall be liable to the end user or consumer for any lack of conformity that exists at the time of delivery of the product under the terms, conditions, and deadlines established by the applicable regulations.

It is recommended that repairs be carried out by qualified personnel.

If at any moment you detect any problem with your product or have any doubt, do not hesitate to contact the official Cecotec Technical Support Service at +34 963 210 728.

9. COPYRIGHT

The intellectual property rights over the texts in this manual belong to CECOTEC INNOVACIONES, S.L. All rights reserved. The contents of this publication may not, in whole or in part, be reproduced, stored in a retrieval system, transmitted, or distributed by any means (electronic, mechanical, photocopying, recording or similar) without the prior authorization of CECOTEC INNOVACIONES, S.L.

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