cecotec

TOOTHCARE JET SPLASH

Irrigador portátil / Portable irrigator

Manual de instrucciones
Instruction manual
Manuel d'instructions
Bedienungsanleitung
Manuale di istruzioni
Manual de instruções
Handleiding
Instrukcja obsługi
Návod k použití
Kullanma kılavuzu

Instrucciones de seguridad	5
Safety instructions	9
Instructions de sécurité	12
Sicherheitshinweise	17
Istruzioni di sicurezza	21
Instruções de segurança	25
Veiligheidsinstructies	29
Instrukcja bezpieczeństwa	33
Bezpečnostní pokyny	38
Güvenlik talimatları	41

ÍNDICE

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SAFETY INSTRUCTIONS

Read these instructions thoroughly before using the appliance. Keep this instruction manual for future reference or new users.

- Ensure that the mains voltage matches the voltage specified on the device rating label and that the plug is earthed.
- This appliance can be used by children aged 8 years and above and people with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children must not play with the appliance. Cleaning and user maintenance should not be carried out by unsupervised children.
- The device must be powered at a low safety voltage as stated on the marking.
- This device should only be used with water or specific solutions.
- To avoid risk of burns, electric shock, fire, or injury:
- Do not attempt to connect or disconnect the charger with wet hands. To disconnect it, pull on the adapter, not the cable.
- Ensure that the mains voltage matches the voltage specified on the device rating label and that the plug is earthed.
- Check the charging cable regularly for visible damage. If the cord or adapter is damaged, malfunctioning, or has been dropped, it must be repaired by Cecotec's official Technical Support Service to avoid any danger.
- Disconnect the charger from the mains when the product is charged.

- Keep the product and cord away from fire and other heat sources.
- Never leave a product unattended while plugged in.
- You must charge the product with the charger supplied by Cecotec. Do not use any other charger type.
- Use the appliance only for the purposes described in this manual.
- Only use tips or accessories recommended or marketed by Cecotec.
- Do not drop or insert any foreign object into any opening or hose.
- Do not direct water under the tongue, ear, nose, or other sensitive areas. This product can produce pressures which may cause serious damage in these areas.
- Do not use the appliance if you are wearing oral jewellery. Remove it before use.
- Do not use the appliance if you have an open wound on your tongue or your mouth. If you have been advised by your physician to receive antibiotic pre-medication prior to dental procedures, you should consult your dentist before using this device or any other oral-hygiene device.
- This appliance is designed for domestic use only and is not intended for bars, restaurants, farmhouses, hotels, motels, and offices.
- Store the appliance and its instruction manual in a safe and dry place when not in use.
- The oral irrigator may cause bleeding gums when first used. This is completely normal and will stop happening with frequent use. Adjust the intensity of the appliance to reduce the presence of bleeding gums.
- Do not use the appliance without assembling a tip.
- Do not place the appliance in a hot or wet place, nor in places with direct sunlight.

- During cleaning or maintenance, make sure that the appliance is switched off and is not charging.
- Users suffering from periodontitis or those who have undergone oral surgery within the last 2 months should consult a physician before using this appliance.
- Do not try to open, disassemble, or repair the appliance by yourself. There are no serviceable parts inside the device.
- To avoid burns, do not fill the tank or attempt to clean your mouth with water above 40 °C.
- Do not use the device for any purpose other than oral hygiene.
- To avoid injury, do not press the tip between your gums or teeth.
- Do not use iodine, salt water, or concentrated essential oils that are insoluble in water with this appliance. Doing so would reduce performance and shorten its shelf life.
- Do not fill the device with pure mouthwash. If you want to use mouthwash, dilute it according to the appropriate proportions.

Battery instructions

- This device contains an irreplaceable battery.
- This device includes a Li-ion battery. Do not incinerate or expose it to high temperatures, as it may explode.
- WARNING: to recharge the battery, use only the power cord supplied with the appliance.
- Leaks from the battery or cells can occur under extreme conditions. Do not touch any liquid that leaks from the battery. If the liquid gets into contact with skin, wash immediately with soap and water. If the liquid gets into the eyes, wash them immediately with clean water for a minimum of 10 minutes and seek medical attention. Wear gloves to handle

the battery and dispose of it immediately in accordance with local regulations.

- Keep batteries out of the reach of children when there is a risk of ingestion.
- Battery ingestion can cause burns, soft tissue perforation, and death. It can cause severe burns within two hours of the ingestion.
- In case of battery ingestion, please seek medical help immediately.
- Do not expose batteries to heat or fire. Avoid storage in direct sunlight.
- Do not subject batteries to mechanical shock.
- Do not leave the battery on prolonged charge when not in use.
- After long periods of storage, it may be necessary to charge and discharge the batteries several times to obtain maximum performance.
- Dispose of the battery properly.
- The recommended room-temperature range during charging is between 4 °C and 40 °C.
- The recommended room temperature range for battery storage is between 0 to 50 °C.
- The recommended room-temperature range during charging is between 5 and 50 °C.

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1. PARTS AND COMPONENTS

Fig. 1

1. Nozzle
2. Water inlet plug
3. Nozzle release button
4. Power button
5. Mode change button
6. Type C charging port
7. Water tank
8. Nozzle
9. Charging cable

NOTE:

The graphics in this manual are schematic representations and may not exactly match the device.

2. BEFORE USE

- This appliance is packaged in a way as to protect it during transport. Take the appliance out of its box and remove all packaging materials. You can keep the original box and other packaging elements in a safe place to prevent damage to the device if you need to transport it in the future. In case the original packaging is disposed of, make sure all packaging materials are recycled accordingly.
- Make sure all parts and components are included and in good conditions. If there is any piece missing or in bad conditions, contact the official Cecotec Technical Support Service immediately.
- Do not remove the product's serial number in order to keep proper traceability if technical assistance is required.

Box content

- ToothCare Jet Splash
- 2 nozzles
- Charging cable
- Travel bag
- Instruction manual

3. OPERATION

Make sure that the battery is fully charged before using the appliance for the first time. You cannot use the appliance while it being charged.

Filling up the water tank

Make sure that the appliance is switched off before adding water to the tank.

Fig. 2

1. Remove the water inlet plug by pulling it upwards and fill the tank.
2. You can also fill the water tank by removing it from the main body. To do this, pull the water tank downwards to remove it, fill the water tank under the tap and replace it.

Fitting the nozzle

Fig. 3

Insert the nozzle into the nozzle hole. You will hear it click into place if correctly installed.

Turn the nozzle to set the water outlet direction.

To remove the nozzle, press the release button and pull it upwards.

Selecting modes

Fig. 4

You can select the mode of use from the 3 available before switching on the appliance. To do this, press the mode change button repeatedly. The indicator light for each mode will light up for 5 seconds.

- The top light corresponds to the gentle mode.
- The central light corresponds to the normal mode.
- The bottom light corresponds to the DIY mode.

Note: the DIY mode pressure can only be adjusted when the irrigator is switched on.

You can also change the mode while using the irrigator by repeatedly pressing the mode change button.

Setting the water pressure in the DIY mode

With the appliance switched on and the DIY mode selected, press and hold the mode change button for 1 second. The DIY mode indicator light will flash. Release the button, the irrigator will change the pressure between the 8 available levels. It will hold each level for 1 second and move on to the next level; when it reaches level 8, it will start the cycle again. When the irrigator is at the desired pressure level, press the mode change button once to select the desired mode, the DIY mode indicator light will return to steady.

The irrigator will save this level as the default for DIY mode until you manually change it again.

Using the irrigator

1. Fill the water tank and insert a nozzle into the irrigator.

ENGLISH

2. Select the desired mode and insert the irrigator nozzle into your mouth.

NOTE: when using the appliance for the first time, start with the Gentle mode or DIY mode at its lowest level. Water pressure can increase gradually as your mouth adapts to it. You may experience itching or slight bleeding of the gums during use. This is normal and it should disappear after 1-2 weeks of continuous use.

3. Point your mouth towards the sink and keep your mouth ajar to let the water drain out.
4. Press the power button once to switch the appliance on.
5. Start cleaning from the back molars and move the irrigator, keeping it at 90°, along the gum line and between the teeth until it reaches the incisors. Continue to the molars on the other side.

NOTE: make a continuous movement, do not stop too long on each tooth so as not to cause damage to the gums.

6. Repeat this step on the outer and inner sides of each jaw (Fig. 5).
7. If you need to refill the water tank during use, switch off the irrigator, refill the tank and then turn it back on to continue cleaning.
8. When you have finished using the appliance, switch it off by pressing the power button.

The appliance will automatically switch off after 2 minutes of continuous operation.

Note: the irrigator memorises the last setting, so when you switch it on again it will automatically programme the last mode used.

9. Empty the remaining water from the tank and dry it after each use.

Note: the irrigator has a protection function which means that when the voltage exceeds 4.5 A, the appliance switches off automatically.

Checking the battery charge level

With the irrigator switched off, press and hold the mode change button to check the battery level. If there is 40% or less battery remaining, the DIY mode indicator light will light up; if there is between 80% and 41% battery remaining, the DIY and Normal mode indicator lights will light up; if the battery is between 100% and 81%, the DIY, Normal and Gentle mode indicator lights will light up.

If during use, the indicator light of the selected mode starts flashing, it means that the battery level is low.

Charging the irrigator

Connect the charging cable to the appliance and then to a compatible power adapter.

While the appliance is charging, the DIY mode indicator light will flash slowly. When fully charged, the 3 mode indicator lights will stay on for approximately 20 minutes and then turn off.

4. CLEANING AND MAINTENANCE

Before cleaning the appliance, make sure that it is switched off.

Fig. 6.

- Remove the water inlet plug. Turn the irrigator upside down to remove all water.
- Remove the water tank by pulling it and rinse its interior with clean water at room temperature. Do not clean the tank with any type of cleaning agent.
- Remove the nozzle and clean it with water.
- Only use clean water or a neutral cleaning agent to clean the body of the appliance. Do not use abrasive or corrosive products. The temperature of the water used to clean the appliance must not exceed 50 °C, otherwise the appliance may get damaged.
- Do not immerse the body of the appliance in water. If necessary, clean the surface of the appliance with a soft cloth and a little neutral soap.
- After cleaning the irrigator, dry all parts thoroughly before storing.
- Store the appliance in a dry place away from sunlight.

Removing limescale deposits

Depending on the mineral content of the water used, limescale deposits may build up in the appliance. Failure to remove them will reduce the effectiveness of the appliance.

- Point the nozzle towards the sink, switch on the appliance and let all the water flow out of the tank.
- Fill the water tank with clean, warm water and repeat the previous step to clean the tank. This cleaning process should be repeated every 1-3 months.

5. TROUBLESHOOTING

Problem	Possible causes	Possible solutions
There is a water leak between the nozzle and the body of the appliance	The nozzle is not correctly fitted.	Fit the nozzle correctly.
There is a leak in the water inlet	The water inlet is dirty.	Clean the water inlet.
	The water tank plug is not correctly positioned.	Fit the plug correctly.
There is a leak in the water tank	The tank is not assembled correctly.	Position the tank correctly.
	The rubber seal is damaged.	Contact the official Cecotec Technical Support Service.

ENGLISH

Water pressure is insufficient	The selected pressure level is low.	Select a higher pressure level.
	Battery level is low.	Charge the appliance.
The appliance does not switch on	Battery level is low.	Charge the appliance and switch it on again.
The appliance does not expel water	It is the first time you use it, or you have not used it for a long time.	Ensure there is water in the tank. Make sure that there are no limescale deposits in the nozzle. Switch on the appliance and let it run for a few seconds until water starts to come out.

6. TECHNICAL SPECIFICATIONS


Product: ToothCare Jet Splash
Product reference: EU01_100061

Product input:
Voltage: 5 V =
Intensity: 1 A
Battery (not replaceable): Lithium-ion, 3.7 V, 1100 mAh, 4.07 Wh
IPX7

To recharge the battery, we recommend using an adapter with the following specifications:
Input: 100-240 V~, 50/60 Hz
Output: 5 V =, 1 A

Technical specifications may change without prior notification to improve product quality.
Made in China | Designed in Spain

7. DISPOSAL OF OLD ELECTRICAL AND ELECTRONIC APPLIANCES

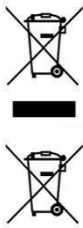This symbol indicates that, according to the applicable regulations, the product and/or batteries must be disposed of separately from household waste. When this product reaches the end of its shelf life, you should dispose of the cells/batteries/accumulators and take them to a collection point designated by the local authorities.

Consumers must contact their local authorities or retailer for information concerning the correct disposal of old appliances and/or their batteries. Compliance with the above guidelines will help protecting the environment.

You can find information about national packaging recycling systems and their labelling on our website.

8. TECHNICAL SUPPORT AND WARRANTY

Cecotec shall be liable to the end user or consumer for any lack of conformity that exists at the time of delivery of the product under the terms, conditions, and deadlines established by the applicable regulations.

It is recommended that repairs be carried out by qualified personnel.

If at any moment you detect any problem with your product or have any doubt, do not hesitate to contact the official Cecotec Technical Support Service at +34 963 210 728.

9. COPYRIGHT

The intellectual property rights over the texts in this manual belong to CECOTEC INNOVACIONES, S.L. All rights reserved. The contents of this publication may not, in whole or in part, be reproduced, stored in a retrieval system, transmitted, or distributed by any means (electronic, mechanical, photocopying, recording or similar) without the prior authorization of CECOTEC INNOVACIONES, S.L.

10. SIMPLIFIED DECLARATION OF CONFORMITY

Cecotec Innovaciones hereby declares that this appliance complies with the essential requirements and other relevant provisions of the regulations applicable in the European Union. This appliance has been designed, manufactured and tested in compliance with required safety and quality standards. The full text of the EU Declaration of Conformity can be found on the following website: <https://cecotec.es/es/information/declaration-of-conformity>

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