
cecotec

TOOTH CARE ULTRASONIC ELITESWING

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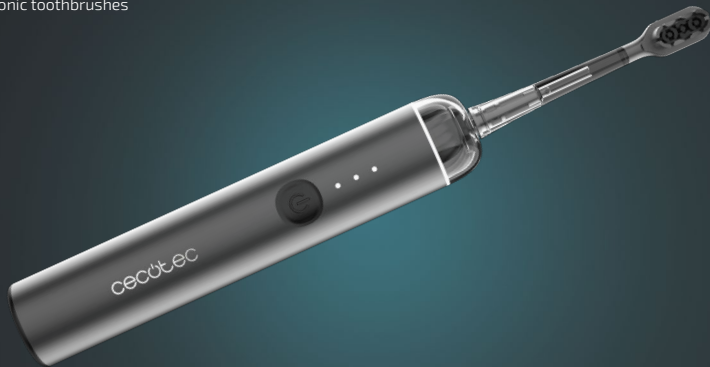SAFETY INSTRUCTIONS

Read the following instructions carefully before using the appliance. Keep this instruction manual for future reference or new users.

- Make sure that the mains voltage matches the voltage specified on the product's rating label and that the plug is earthed.
- This appliance can be used by children aged 8 years and above and people with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children must not play with the appliance. Cleaning and user maintenance must not be carried out by

children without supervision.

- The appliance must be powered at a low safety voltage as stated on the marking.
- This device should only be used with water or specific solutions.
- To avoid risk of burns, electric shock, fire, or injury:
- Do not attempt to connect or disconnect the charger with wet hands. To disconnect it, pull on the adapter, not the cable.
- Ensure the mains voltage matches the voltage specified on the appliance's rating label and that the socket is properly earthed.
- Check the charging cable regularly for visible damage. If the cord or adapter is damaged, malfunctioning, or has been dropped, it must be

repaired by Cecotec's official Technical Support Service to avoid any danger.

- Disconnect the charger from the mains when the product is charged.
- Keep the product and cord away from fire and other heat sources.
- Never leave a product unattended while plugged in.
- You must charge the product with the charger supplied by Cecotec. Do not use any other charger type.
- Use the appliance only for the purposes described in this manual.
- Only use tips or accessories recommended or marketed by Cecotec.
- Do not drop or insert any foreign object into any

opening or hose.

- Do not direct water under the tongue, ear, nose, or other sensitive areas. This product can produce pressures which may cause serious damage in these areas.
- Do not use the device if you are wearing oral jewellery. Remove it before use.
- Do not use the device if you have an open wound on your tongue or your mouth. If you have been advised by your physician to receive antibiotic pre-medication prior to dental procedures, you should consult your dentist before using this device or any other oral-hygiene device.
- This appliance is designed for domestic use only and is not intended for bars, restaurants, farmhouses, hotels, motels, and offices.

- Store the appliance and its instruction manual in a safe and dry place when not in use.
- The oral irrigator may cause bleeding gums when first used. This is completely normal and will stop happening with frequent use. Adjust the intensity of the appliance to reduce the presence of bleeding gums.
- Do not use the device without assembling a tip.
- Do not place the device in a hot or wet place, nor in places with direct sunlight.
- During cleaning or maintenance, make sure that the device is switched off and is not charging.
- Users suffering from periodontitis or those who have undergone oral surgery within the last 2 months should consult a physician before using

this device.

- Do not try to open, disassemble, or repair the appliance by yourself. There are no serviceable parts inside the appliance.
- To avoid burns, do not fill the tank or attempt to clean your mouth with water above 40 °C.
- Do not use the device for any purpose other than oral hygiene.
- To avoid injury, do not press the tip between your gums or teeth.
- Do not use iodine, salt water, or concentrated essential oils that are insoluble in water with this appliance. Doing so would reduce performance and shorten its shelf life.
- Do not fill the device with pure mouthwash

liquid. If you want to use mouthwash, dilute it according to the appropriate proportions.

Battery instructions

- This appliance contains an irreplaceable battery.
- This appliance includes a Li-ion battery. Do not incinerate or expose it to high temperatures, as it may explode.
- WARNING: to recharge the battery, use only the power cord supplied with the appliance.
- Leaks from the battery or cells can occur under extreme conditions. Do not touch any liquid that leaks from the battery. If the liquid comes into contact with skin, wash immediately with soap and water. If the liquid gets into the eyes,

wash them immediately with clean water for a minimum of 10 minutes and seek medical attention. Wear gloves to handle the battery and dispose of it immediately in accordance with local regulations.

- Keep batteries out of the reach of children when there is a risk of ingestion.
- Battery ingestion can cause burns, soft tissue perforation and death. Severe burns may occur within two hours of ingestion.
- In case of battery ingestion, please seek medical help immediately.
- Do not expose batteries to heat or fire. Avoid storage in direct sunlight.
- Do not subject batteries to mechanical shock.
- Do not leave the battery charging for an

extended period of time when not in use.

- After long periods of storage, it may be necessary to charge and discharge the batteries several times to obtain maximum performance.
- Dispose of the battery properly.
- The recommended room-temperature range during charging is between 4 °C and 40 °C.
- The recommended room temperature range for battery storage is between 0 to 50°C.
- The recommended room-temperature range during charging is between 5 and 50 °C.

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1. PARTS AND COMPONENTS

Fig. 1

- 1. Brush head
- 2. Warning light
- 3. Gum Care mode
- 4. Cleaning mode
- 5. Whitening mode
- 6. High-frequency mode
- 7. Power button
- 8. Handle
- 9. Charging port

NOTE:

The graphics in this manual are schematic representations and may not exactly match the product.

2. BEFORE USE

- This appliance comes in a packaging designed to protect it during transport. Take the appliance out of its box and remove all packaging materials. You can keep the original box and other packaging materials in a safe place to prevent damage to the appliance if you need to transport it in the future. If you wish to dispose of the original packaging, make sure all items are recycled properly.
- Check that all parts and components are included and in good condition. If any of them are missing or damaged, please contact Cecotec's Official Technical Support Service immediately.

Box content

- ToothCare UltraSonic EliteSwing
- 2 brush heads
- Charging cable
- Instruction manual

3. OPERATION

- The electric toothbrush has undergone a complete discharge and charge cycle before leaving the factory. However, due to inventory time variations and normal natural discharge, the battery in your newly purchased product may not be fully charged. It is recommended to recharge it before use.
- Ensure your local power supply voltage matches the voltage range indicated on the product's rating plate. Connect the electric toothbrush to the charging cable, select a suitable charging head, and plug it into a secure power outlet.
- When the battery is low, a red light flashes at the top to indicate that the electric toothbrush needs recharging.
- Charging takes approximately 5 hours, and the battery level LED indicator light flashes during charging.
- Once fully charged, the battery level LED indicator light stops flashing and remains steadily lit.
- After being fully charged, promptly remove the charging head from

the socket; do not leave it continuously plugged in.

Brush head installation

When installing the brush head, align the bristle front end with the handle front. Align and insert the drive shaft into the hole at the bottom of the brush head, push firmly until fully seated, leaving a 0.2-1 mm gap between the brush head bottom and handle top.

Fig. 2

Functions

- 1. Timer function: In each operating mode, the brush head pauses briefly (0.5 seconds) every 30 seconds to prompt users to change brushing position, automatically stopping after 4 sets of 30 seconds (approximately 120 seconds total).
- 2. Pressure sensor function: During brushing, the top indicator light shows green under normal pressure. If excessive pressure is applied, the top indicator light turns red to remind users to adjust the pressure.

- 3. Auto shut-off function: When the device remains inactive for 10 seconds in standby mode, it automatically enters auto shut-off mode and the indicator light turns off. As soon as the user picks up or shakes the toothbrush, it automatically activates, the indicator light illuminates, and the brush head position calibrates.
- 4. Activation function: When the user picks up or shakes the toothbrush, it automatically activates, the indicator light illuminates, and the brush head position calibrates.
- 5. Travel lock function: When the device is switched on, press and hold the button for more than 5 seconds to activate the travel lock. The motor will sound an alert and the light panel will display a red light. In lock mode, pressing the button once causes both the mode light and top red light to flash for 5 seconds before turning off. To unlock, press and hold the button for over 3 seconds or press twice. The motor sounds when unlocking.
- 6. Charging status: When battery is low, the third indicator light slowly flashes. At medium battery level, the third and second lights cycle up and down. When nearly full, all three indicator lights cycle up and

down; when fully charged, all three indicator lights remain steady. The device cannot be switched on while charging.

- 7. Low voltage protection function: When voltage falls below 3.5 V, the top panel's red light flashes for 5 seconds indicating insufficient power; when voltage falls below 3.3 V (to protect against over-discharge), the device automatically prevents operation.
- 8. Incorrect discharge prevention function: When in standby mode, if continuously shaken for over 30 seconds, the device automatically enters auto shut-off mode and turns off all indicator lights. To resume, press the button once to return to standby mode. (This prevents battery drainage if users forget to activate the travel lock during transport).

Mode selection

- 1. Gum Care mode: Top light on
- 2. Cleaning mode: Middle light on
- 3. Whitening mode: Bottom light on
- 4. High-frequency mode: All three lights on

Using the toothbrush

- Moisten the brush head and apply a small amount of toothpaste.
- Place the brush head in your mouth with the bristles at a slight angle (45 degrees) against your teeth. Press the power button to activate brushing mode and the toothbrush will begin operating. Continuously press the power button to cycle through four modes. When switching to each brushing mode, the corresponding indicator light turns white (as shown in Figure 3).
- Gently press the brush head against your teeth and brush each quadrant for 30 seconds. If brushing stops automatically, simply press the power button to turn it off. Every 30 seconds, there will be a brief pause to remind you to change cleaning area (as shown in Figure 4).
- After completing all four quadrants, the brushing mode automatically stops. Remove the brush head from your mouth to

prevent toothpaste foam splatter. (Within 5 seconds of switching to brushing mode, press and hold the power button to stop brushing mode).

Fig. 3 Fig. 4

4. CLEANING AND MAINTENANCE

- User cleaning and maintenance tasks must not be performed by children. Before rinsing with clean water, please turn off the device and remove the brush head. It is recommended to replace the brush head every 3 months.
- Rinse the brush head with clean water, shake off excess water after cleaning, and do not use water hotter than 80°C for cleaning.
- Rinse the body with clean water, do not scrub or damage the waterproof rubber seal. Do not immerse this device in water for extended periods. Do not use water hotter than 40°C for cleaning.
- Dry with a clean cloth and store in a dry, safe place.

Fig. 5.

5. TROUBLESHOOTING

Problem	Possible cause	Solution
The brush head does not move	Low battery	Charge the toothbrush
The brush head vibrates weakly		
Short operating time until it stops working	The battery has reached the end of its useful life	Charge the toothbrush
	Insufficient charging period	Charge the device for at least 5 hours

6. DISPOSAL OF OLD ELECTRICAL AND ELECTRONIC APPLIANCES

This symbol indicates that, according to the applicable regulations, the product and/or batteries must be disposed of separately from household waste. When this product reaches the end of its service life, you must remove the batteries/accumulators and take it to a collection point designated by local authorities.

Consumers must contact their local authorities or retailer for information concerning the correct disposal of old appliances and/or their batteries. Compliance with the above guidelines will help protecting the environment.

7. BATTERY RECYCLING

- Remove the bottom cover of the product and the two exposed screws.
- Using pliers, grip the upper part of the shaft and pull upwards.

- Remove the battery from its holder and cut the wires connecting it to the PCB.

8. COPYRIGHT

The intellectual property rights over the texts in this manual belong to CECOTEC INNOVACIONES, S.L. All rights reserved. The content of this publication may not, either in part or in its entirety, be reproduced, stored in a retrieval system, transmitted or distributed by any means (electronic, mechanical, photocopying, recording or similar) without prior authorisation from CECOTEC INNOVACIONES, S.L.

9. SIMPLIFIED EU DECLARATION OF CONFORMITY

Cecotec Innovaciones hereby declares that this product complies with the essential requirements and other relevant

provisions of the regulations applicable in the European Union. This appliance has been designed, manufactured and tested in compliance with required safety and quality standards. The full text of the EU Declaration of Conformity can be found on the following website: <https://cecotec.es/es/information/declaration-of-conformity>

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- Img. 1
- 1. Brosse
 - 2. Voyant d'avertissement
 - 3. Mode Gum care
 - 4. Mode Cleaning
 - 5. Mode Whitening
 - 6. Mode High frequency
 - 7. Bouton marche/arrêt
 - 8. Poignée
 - 9. Port de charge

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